

REJUUVENA™

7000

**REFRESH AND REJUVENATE
A SPA LIKE
TREATMENT AT
HOME!**

Look younger by smoothing out and revitalizing your appearance with the Rejuvena 7000.

Experience it in 2 steps - Gentle Soft Peeling and Ion Treatment with Red Light Heat Lamp.

When our skin is unable to exfoliate naturally, it accumulates dirty old cells and its appearance becomes dull. Richway & Fuji Bio has scientifically developed the Rejuvena 7000 and created it to respond to the real needs of aging skin. Let the Rejuvena 7000 gently and naturally reveal your smooth lustrous skin.

**REJUUVENA™ 7000
HI-TECH ION TREATMENT
FACIAL CARE**



THE
RICHWAY
GROUP

RICHWAY & FUJI BIO INC.

1314 South King Street Suite 520, Honolulu, HI 96814 USA Tel: 808-589-2800 Toll-Free: 855-338-6410 Fax: 808-597-1651


www.richwayandfujibio.com www.richwaykorea.com

SHAOXING RICHWAY TRADING CO., LTD., RICHWAY EUROPE B.V., RICHWAY & LIFE JAPAN CORPORATION, R & L CORPORATION

SOFT PEELING MODE

The Soft Peeling mode is painless and gentle on your skin. It uses 30,000 vibrations per second to deeply cleanse and remove contaminants and old cells, revealing your natural skin.

✓ HOW TO USE THE SOFT PEELING MODE

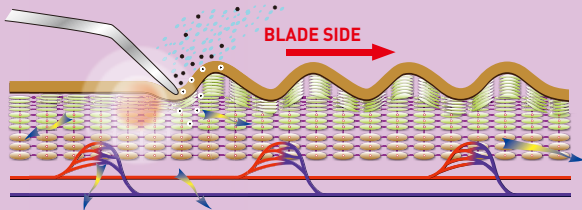

- 1 Wash your face and remove any make up. For best results, steam your skin with a hot towel or facial steamer.
- 2 Use the supplied bottle and mist water on your face. For best results, we recommend using water with a pH of 6.5.
- 3 Turn the unit on by pressing the Power button once. Press it once more for the Soft Peeling mode.
- 4 Gently apply the blade to your skin (blade side), according to directions in Step 1 (Soft Peeling Function).
- 5 By slowly moving the blade against the skin, it removes old cells and other toxic substances (Use for 5 minutes).
- 6 When you are finished, dab your face with a cold wet towel.
- 7 Turn the unit off. (Press the Power button once to change it to Ion Treatment mode. Press once more and hold for 3 seconds to turn off.)

► **NOTE:** After 5 minutes, it will automatically turn to Standby mode. If you want to keep using it on Soft Peeling mode, press the button 1 more time.


Caution (Soft Peeling Mode)



- Use twice a week for normal skin, once in a week for sensitive skin.
- Do not use the Soft Peeling mode every day.
- You may feel a slight tingling after 5 minutes of use.

REJUVENA™ 7000

Combining Soft Peeling massage to enhance the effectiveness of advanced technology.

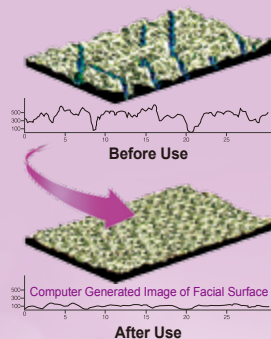

STEP 1 SOFT PEELING FUNCTION [BLADE SIDE]

CLEANSE, TONE, TIGHTEN, AND
LIFT AGING SKIN.

FACIAL CARE SPRAY

We included an empty bottle in the Rejuvena 7000 kit. Fill this with water and apply it to your skin before using the Soft Peeling mode. We recommend water with a pH of 6.5 for best results.

- Use the soft peeling mode two times a week. For users with sensitive skin, use once a week. Use a hot towel or facial steamer to help remove old cells.



1 Facial Care Spray

Prepare your skin by applying water to it. Re-apply water as needed.



2 Start by Gently Massaging Forehead and Cheeks

Gently stroke the forehead and cheeks, follow the directional arrows 3 times on each side.



3 T-zone & Sides of the Nose

Stroke the T-zone & sides of the nose in an upward motion 10 times.



4 Mouth & Cheek Area

Stroke the area around the mouth in a forward motion 10 times.



5 Lines Around Cheek, Nose, Above & Under Lips

Stroke in an upward motion 10 times.



6 Neck Area

Gently stroke lines on neck in an upward motion 10 times. (Do not touch center area)

These statements have not been evaluated by the Food and Drug Administration. These products are not intended to diagnose, treat, cure or prevent any disease. Specific medical advice should be obtained from a licensed health care practitioner.

STEP 2

ION TREATMENT & RED LIGHT HEAT LAMP FUNCTION [FLAT SIDE]

PURE C SERUM

Our advanced Pure C Serum is a unique, oil-free patented Vitamin C formulation that contains powerful antioxidants designed to help reduce visible signs of aging and sagging skin. Apply before using the Rejuvena in Ion Treatment mode.

Available for separate purchase ▶

▶ Use the Ion Treatment mode every day or even twice a day.



1 Ionized Vitamin C

Apply Vitamin C Serum to your skin for a healthy, youthful glow. It will help to reduce the appearance of fine, dry skin lines. (Item not included, must order separately)



2 Start by Gently Massaging Forehead & Cheeks

Stroke the forehead and cheeks in an outward motion 5 times. Avoid the eye area.



3 T-zone & Sides of the Nose

Stroke the T-zone & sides of the nose in an upward motion 10 times.



4 Mouth & Cheek Area

Stroke the area around the mouth in an outward motion 10 times.



5 Lines Around Cheek, Nose, Above & Under Lips

Stroke in an upward motion 10 times.



6 Neck Area


Gently stroke lines on neck in an upward motion 10 times. (Do not touch center area)

These statements have not been evaluated by the Food and Drug Administration. These products are not intended to diagnose, treat, cure or prevent any disease. Specific medical advice should be obtained from a licensed health care practitioner.

ION TREATMENT MODE

The Ion Treatment mode uses a current to increase the efficacy of skin nourishing products by deeply pushing nutrition into your skin.

✓ HOW TO USE THE ION TREATMENT MODE WITH RED LIGHT HEAT LAMP

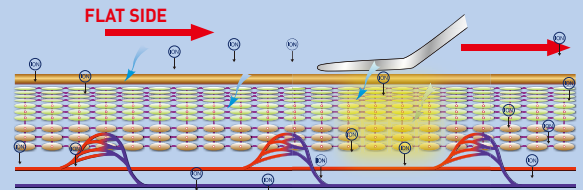

- 1 You can use the Ion Treatment mode with or without using the Soft Peeling mode.
- 2 After the Soft Peeling massage, prepare your skin with Vitamin C Serum. If you do not have Vitamin C Serum use a water-based skin ointment.
- 3 Turn on the unit by pressing the Power button once. Press it twice more for the Ion Treatment mode with Red Light Heat Lamp.
- 4 Gently apply the blade to your skin (flat side), according to directions in Step 2 (Ion Treatment & Red Light Heat Lamp Function).
- 5 By slowly moving the blade against the skin, it helps nutrients penetrate deeper into the skin. (Use for 5 minutes)
- 6 When you are finished, apply an essential cream.
- 7 Turn the unit off. (Press the Power button once to change it to Standby mode. Press once more and hold for 3 seconds to turn off.)

▶ **NOTE:** After 5 minutes, it will automatically turn to Standby mode. If you want to keep using it on Ion Treatment mode, press the button 2 more times.

Caution

(Ion Treatment Mode)




- Be sure to apply Vitamin C Serum or water-based skin ointment before using Ion Treatment mode.
- Use regularly if your skin is dry.
- To avoid excessive use, the device powers off every 5 minutes.
- If you need to, you can select the mode and use it again.

660 nm Red Light Heat Lamp

When using Ion Treatment mode, the Red Light Heat Lamp aids in the absorption of nutrients into the skin.

Red Light

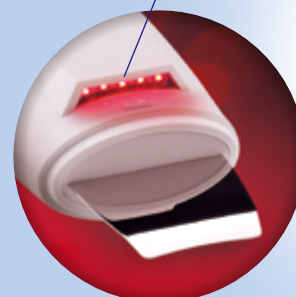

WHERE WONDERFUL SKIN BEGINS!

REJUVENA™ 7000



• 30,000 vibrations per second

Blade Side (Soft Peeling Mode)

Use it with water to remove waste products such as excess skin keratin, sebum, etc.

Flat Side (Ion Treatment Mode)

Use it with cosmetics to increase absorption of active ingredients. The red light heat lamp on the other side helps cosmetics to absorb better during Ion Treatment mode.

Ion Treatment with Red Light Heat Lamp

Soft Peeling Lamp

Power/Charging Lamp

Power Button

Hold the button for 3 seconds. It will make a beeping sound and turn off.

Hand Electrode

Charging Cable Connect

PRODUCT STANDARD

- ▶ **MODEL:** REJUVENA 7000
- ▶ **VOLTAGE:** DC 5V
- ▶ **BATTERY:** 3.7V / 450mAh Li-polymer
- ▶ **MATERIAL:** ABS / PC / etc.

- ▶ **CHARGING TIME:** About 90 min
- ▶ **USE TIME:** About 60 min
- ▶ **SIZE:** 169 x 46 x 25 mm
- ▶ **WEIGHT:** About 95g

▶ INSTALLATION


1. First, find a good location to place the Rejuvena 7000. (A flat and stable surface)
2. Insert the connector of the charger into the Rejuvena 7000.
3. Plug the charger into an electrical socket. When the unit is connected properly, the indicator will light, and charging will begin.

▶ HOW TO USE (Mode Selection)

1. Press Power button once to turn it on (Standby Mode).
2. Press once more, to turn on "Soft Peeling Mode".
3. Press once more, to turn on "Ion Treatment Mode with Red Light Heat Lamp".
4. Press and hold for 3 seconds to turn off.

▶ HOW TO CHARGE

1. You can charge the unit by plugging in the USB cable (included) or a charging cable for a Smartphone (MHL cable).
2. Power/Charging Lamp will flash on & off while charging.
3. Once fully charged, the lamp will stop flashing.
4. Recommended to use over 5V/1A capacity items when using an adapter.

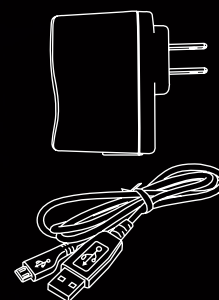

▶ HOW TO CLEAN/STORE?

Use a soft, clean cloth or tissue paper to wipe the blade after each use. Alcohol can be used to disinfect the blade. Avoid letting water get inside of the machine during use and be sure to remove any water after use. Do not place it in water or use soap to clean. Avoid hot or moist places and store it in a cool, dry place.

BE CAREFUL WHEN USING THE REJUVENA™ 7000

- ▶ Keep it away from water in the bathroom or sauna.
- ▶ When used for more than 5 minutes, the unit will turn off automatically. This is a safety feature.
- ▶ If your skin becomes irritated, stop using immediately and consult your therapist.
- ▶ Do not use on irritated skin, such as abscesses or eczema.
- ▶ When finished massaging, clean the blade with a soft, clean cloth or tissue paper.
- ▶ The Rejuvena 7000 is a high-tech precision instrument. Do not store it where the temperature or humidity is high. Prevent physical shock such as dropping. Do not leave in your car on a hot day.
- ▶ If you have not used the Rejuvena 7000 for an extended period of time, the unit should be charged. It can be used for a total of approximately 2 hours before it needs to be charged again.
- ▶ Keep the charger clean by wiping it with a clean cloth or tissue paper.
- ▶ Only use an adapter that is compatible with the Rejuvena 7000, or the unit can be damaged.
- ▶ If you encounter any problems while using the Rejuvena 7000, contact Richway & Fuji Bio at (808) 589-2800 or 855-338-6410 (toll-free), and ask for the Service Department.
- ▶ Damage caused by mishandling or carelessness will not be covered by the warranty.
- ▶ The warranty will not be honored if you disassemble the instrument.

These statements have not been evaluated by the Food and Drug Administration. These products are not intended to diagnose, treat, cure or prevent any disease. Specific medical advice should be obtained from a licensed health care practitioner.